**ENG** Mediblink gallium thermometer M370

**DESCRIPTION OF THERMOMETER COMPONENTS**

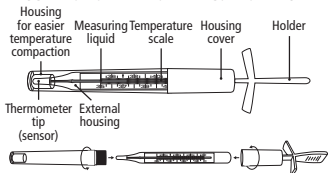
Image 1: Thermometer components

**SPECIFICATION**

Measuring range:	35.5°C to 42.0°C
Resolution:	0.1°C
Accuracy:	±0.1°C/±0.15°C
Operating conditions:	Temperature: -10°C to 42°C (Allows only measuring of temperature above room temperature)
Storage conditions:	Temperature: -15°C to 42°C
Expected use:	Medical glass thermometer with liquid for measuring body temperature.
Description:	Production and testing of the thermometer shall be carried out in accordance with standards.
Quality guarantee:	<ul style="list-style-type: none"> <li>● Directive 93/42/EEC concerning medical devices</li> <li>● Standards EN 12470-1:2000 + A1:2009 - Clinical thermometers, 1st part: Glass thermometers with liquid metal and the display of the highest temperature.</li> </ul>

**⚠ PRECAUTIONS**

1. Make sure the thermometer does not fall or be subjected to sudden impacts.
2. Do not bend or chew the tip of the thermometer.
3. Do not use thermometers with damaged external housing (danger of injury).
4. NOTE: Children should be monitored by an adult throughout the measurement procedure. Keep out of the reach of children.
5. The measuring fluid (a mixture of gallium, tin and indium; hereinafter: "fluid") must not come into contact with precious or light metals as it reacts to them.
6. Care should be taken to avoid any potential risk or damage during measuring.
7. This thermometer cannot be a substitute for a visit to a doctor. If high fever is not reduced, especially in children, the patient should be checked by a doctor.

8. The thermometer is not intended for premature children.

**INSTRUCTIONS ON USE**

1. Prior to measuring the temperature, make sure that the gallium value is below 36°C on the measuring scale.
2. The thermometer shall be placed depending on the method of measuring the temperature (axillary, rectal or oral).
3. The thermometer measures the body temperature for about 5 minutes. You can then read the measured value. Gently turn the thermometer to make it easier to read the temperature value.
4. Reset the liquid below 36°C after measuring has been carried out by returning the thermometer to the housing and making sure the tip of the thermometer is attached to the housing cover (shown in Image 1). Make sure the housing is well closed. Then grab the holder of the housing and shake it with fast rhythmic movements (Image 2) so the liquid falls under the first numbered line.

Image 2: Resetting the temperature value of the thermometer

5. Clean the glass thermometer after use. Please follow the Cleaning instructions carefully.

**⚠ Note:** Always disinfect the thermometer before and after use. Professional staff must always provide for disinfection in hospitals.

**RELIABLE WAYS OF TEMPERATURE MEASURING**

**ORAL** - Measuring time: At least 5 minutes, normal temperature at 37°C. Put the end of the thermometer tip under the tongue and breathe through your nose so the measurements is accurate. Oral measuring is permitted for children aged 10 years or more. Do not chew, talk or play when measuring the temperature in your mouth. Avoid bathing, feeding and smoking for at least half an hour before temperature measuring.

**RECTAL** - Measuring time: At least 5 minutes, normal temperature at 37.5°C. Gently insert the thermometer tip not more than 0.5cm deep into the anal orifice. Stop immediately if you feel resistance. Do not use lubricant. If measuring body temperature of a baby, place the baby on a hard surface and lift its legs almost perpendicularly and hold it by the ankles to keep it still. Insert the thermometer tip very gently into the anal orifice.

**AXILLARY** - Measuring time: Minimum 5 minutes and a normal temperature at 36.7°C. Make sure the armpit is dry. Then put the thermometer tip under the armpit and press your hand against the chest.

**CLEANING INSTRUCTIONS**

The thermometer has hermetically sealed hygienic glass housing and can be disinfected without restrictions. For disinfection, Mediblink M140 alcohol wipes are recommended.

 The temperature of disinfectant and cleaning agents shall not exceed 42°C. Clean the glass thermometer with a cold soap before and after use. Disinfect the thermometer with a cotton pad or cloth soaked in alcohol or immersed in a disinfectant solution (immersion time must be in accordance with the manufacturer's instructions). If the thermometer is shattered, spilled liquid may be removed by a brush, or a cloth, soaked in soapy water. The spilled fluid is non-toxic. Broken parts may be disposed as household waste.

## WARRANTY

### Product: Mediblink gallium thermometer M370

**Manufacturer:** Wuxi Medical Instrument Factory Co., Ltd, No.43 Xixin Road, Zhangjing Xibei Town, Wuxi City, 214194 Jiangsu, China

**Sellers name, address, signature and stamp:**

**Date of extradition/sales:**

## WARRANTY TERMS

Dear customers!

The warranty period is 10 years and starts from the day of product purchase. In case of product claim, you have to show the invoice. We kindly ask you to save the invoice!

Unfortunately, wrong handling with the device is a reason for 95% of customer complains. You can easily avoid any problem, by getting useful information provided by our special service department. To reach our service department, you can call or send e-mail to Mediblink local distributor.

Before sending the product back to retailer, we kindly ask you to call our service department, to get help about how to use the device to save you with unneeded trips.

The manufacturer guarantees free elimination of all imperfections due to defects in material or manufacturing procedure by repairing or replacing the product. In case that the product can not be repaired or replaced, the customer will get the money refund. The guarantee is not valid in case of the force majeure, accidents or unexpected events (such as lightning, water, fire etc.), incorrect use or incorrect transport, non-compliance with safety and maintaining regulations or in case of unprofessional product intervention.

Traces of every day product usage (scratches, abrasions) and not subject to claim. The warranty does not eliminate the customer rights, which originate from seller responsibility for product flaws. By accepting the claimed product by the service department, the service department does not take responsibility for loss of saved data or settings on the product. All product repairs, which are performed out of product warranty period, have to be paid by customer by prior notice.

The manufacturer guarantees the product quality and flawless product operation in the warranty period, which starts with the

day of product purchase. If the product can not be repaired in 45 days, the product will be replaced with a new one. In case that the product can not be replaced, the money will be refunded to the customer.

In case of product claim, call or send e-mail to Mediblink local distributor.

Any serious incident that has occurred in relation to the device should be reported to the manufacturer (Wuxi Medical Instrument Factory Co., Ltd and company Mediblink d.o.o.) and the competent authority of the Member State in which the user and/or patient is established.

## LIST OF SYMBOLS

	Follow the instructions e-IFU: <a href="http://www.mediblink.com/IF/M370.pdf">http://www.mediblink.com/IF/M370.pdf</a>
	Attention
	Group code (mm/yyyy; month/year) marked on the product packaging (written by laser without a LOT series symbol)
	Class Im medical device
	Manufacturer
	Paper
	Number of products in one packaging
	The product does not contain mercury
	Maximum and minimum temperature: Store at a temperature between -10°C and +42°C
	Store in a dry place
	Fragile: handle carefully

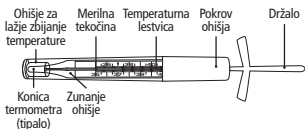
 Wuxi Medical Instrument Factory Co., Ltd  
No.43 Xixin Road, Zhangjing Xibei Town,  
Wuxi City, 214194 Jiangsu, China  
 Lotus NL B.V., Koningin Julianaplein 10,  
1e Verd, 2595AA, The Hague, Netherlands

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